



21.3-inch 3MP Diagnostic Monitor

| 21HQ513D |

Accurate Images for a Precise Review


3MP (1536 x 2048) IPS Display | 1100cd/m² Brightness | Multi-Resolution Mode (2MP/3MP)

21HQ513D has a 21.3-inch 3MP IPS display with 1100cd/m² (Typ.) brightness which facilitates the precise review of CT and Angiography images by enabling medical professionals to detect subtle details faster. High brightness and contrast enhance the visibility of text and help to illuminate abnormalities. The 2MP and 3MP resolution options provide optimized image quality for various modalities, making it easier to select the optimal viewing settings.

Reliable Imagery for Optimized Diagnostics


Focus View | Pathology Mode

Focus view allows users to review specific parts of an image more closely. The area of concern can be selected, magnified, and the brightness can be adjusted depending on the image source, using just a mouse and keyboard, while other areas of the image are automatically darkened. Additionally, with Pathology Mode the 21HQ513D helps medical professionals to see more details with accurate color similar to that shown directly under a microscope.

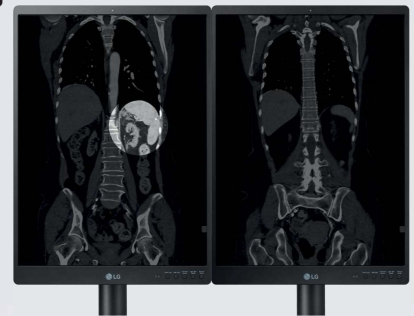

Internal Sensors (Backlight Sensor | Front Sensor | Auto Luminance Sensor)

The 21HQ513D is equipped with a variety of sensors to create an optimized diagnostic environment. The backlight sensor automatically adjusts luminance to stabilize brightness for consistent imagery. The front sensor enables the device to conduct automatic hardware calibration without the need for additional equipment to maintain accurate values and improve image quality and consistency. Not only does this increase efficiency it also reduces the need for additional operating costs or manpower. Finally, the Auto Luminance Sensor ensures screen brightness is always optimized for the ambient lighting conditions, reducing the risk of eye strain.

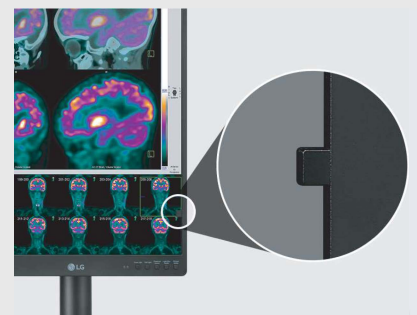

Intuitive Controls for Increased Productivity


5 Hot Keys | Presence Sensor

5 Hot Keys built into the 21HQ513D allow you to access Down/Wall Light, Presence Sensor, Light Box mode, and Picture Mode faster and more intuitively than an on-screen menu. The built-in Presence Sensor automatically turns off the display when no motion is detected. This reduces power consumption while dramatically reducing the hassle of manually turning the display on and off.

Ergonomic Stand (One-click) | Lightweight Body | Down/Wall Light | Narrow Side Bezel

An ultra-lightweight body and one-click stand make installing the 21HQ513D simple. Additionally, users can easily adjust the height, tilt, pivot of the display thanks to its ergonomic design, helping to reduce chronic pain caused by long hours of work and poor posture. Down and Wall Lighting modes reduce the contrast between the display and the ambient surroundings allowing users to work comfortably without having to adjust the lighting to view paper documents in dark reading rooms. Additionally, the Narrow Side Bezel minimizes disruption that may be caused by a wide bezel, meaning that the monitor can be pivoted to create an ideal multi-display setup.

Specification

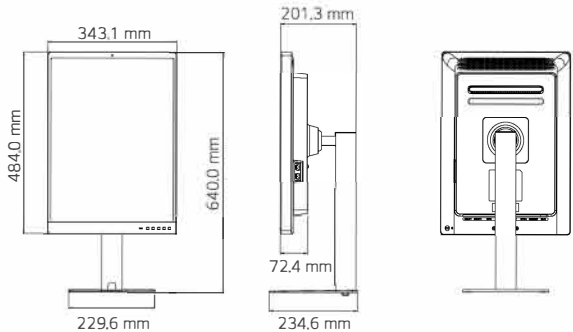Model		21HQ513D
Model		
Panel	Panel Type	IPS
	Inch (Aspect Ratio)	21.3-inch (3:4)
	Resolution	3MP (1536 x 2048)
	Surface Treatment	Anti-glare, 3H
	Color Gamut (Typ.)	NTSC 72%
	Viewing Angles (CR≥10)	178° (Right/Left), 178° (Up/Down)
	Brightness (Typ.)	1100 cd/m ²
	Contrast Ratio (Typ.)	1800:1
	Response Time (GtG*)	5ms
Feature	DICOM Compliant	Yes
	HW Calibration	Yes (PerfectLum) with Internal Front Sensor
	Display Mode	Multi-resolution (3/2MP), Focus View Mode, Light Box Mode, Pathology Mode
	Super Resolution+	Yes
Video Signals	Input Terminals	DisplayPort x2, DVI x1
	Output Terminals	DisplayPort x1 (Daisy Chain, 3MP/ @60Hz)
Connectivity	USB	2 Upstream, 2 Downstream
Power	AC Input	100-240Vac, 50/60Hz
	Power Consumption (Max.)	120W
	Power Consumption (DC Off)	Less than 0.3W
Certifications & Standards		IEC (IEC 60601-1 / IEC 60601-1-2), EN (EN 60601-1 / EN 60601-1-2), IEC (IEC 60950-1 / IEC 55032,55024) EN (EN 60950-1 / EN 55032, 55024), cUL (ANSI/AAMI ES 60601-1, CSA CAN/CSA-C22.2 NO. 60601-1) FCC (FCC part 15 Class A), FDA (510(k)(Class II)), RoHS, REACH, WEEE
User Convenience	Dual Controller	Yes
	Text Mode / Flicker Safe	Yes / Yes
	Lighting	Yes (Down / Wall)
	Hot Key	Yes (5 Keys - Down Light / Wall Light / Presence Sensor / Light Box Mode / Picture Mode)
	Presence Sensor	Yes
	Auto Luminance Sensor	Yes
Physical Specification	Weight (with stand)	7.0 kg (15.4 lbs)
	Weight (without stand)	5.0 kg (11.0 lbs)
	Adjustable Stand	Tilt : -5°~20°, Pivot : 90°, Height Range : 110 mm

*GtG : Gray to gray response time.

**Specifications may vary by region and all features, standards, and other product specifications are subject to change without notice or obligation.



LG Electronics Inc.

<https://www.lg.com/global/business/virtual-showroom/medical-display>

Copyright © 2021 LG Electronics. All Rights Reserved.